BEAUTIFLY

B-Looky

PRO
series

Fale Radiowe

Rejuvenation Radio Frequency

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User manual

/beautifully.eu

@beautifully.eu

/beautifulyeu

beautifully.eu

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RADIO WAVES	STRONG ANTI-AGING ACTION
RED PHOTO LIGHT	12 000 RPM VIBRATION

IT SMOOTHES AND FIRMS THE SKIN AROUND THE EYES

IT HAS A POSITIVE EFFECT ON EDEMA

B—Looky PRO is a powerful anti-aging device for the care of the skin around the eyes. Effectively soothes tired and swollen eyes. It stimulates the skin around it, firms it and smoothes it. It has the function of radio waves with a strong lifting effect. High frequency 12 000 RPM vibrations combined with photon red light perfectly smooth and care for the thin and sensitive skin around the eyes.

Technical data B—Looky PRO

Input:	DC5V/1A
Battery:	500mAh
Unit weight:	54g
Size:	149.9x28 x22mm
Vibration frequency:	12000RPM
RF:	400KHZ
Charging time:	3 h
Using time:	>30min



Part description

Device

1. Massage head
2. Red light
3. LCD display
4. Level button
5. Power / Mode button
6. USB charging port

LCD display

1. Battery indicator
2. LED light indicator
3. Radio frequency indicator
4. Vibration indicator
5. RF 3 working level

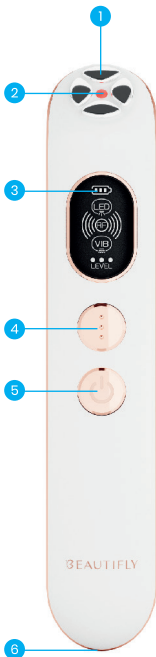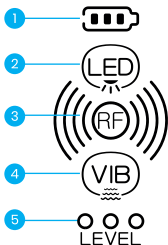Functions

1. B-Looky PRO is a device for the care of the eye area. It uses the frequency of radio waves along with vibrations, and additionally has the function of red photon light. Use it with your eye cream or gel to make your skin smoother, firmer and brighter.
2. Anti-aging technology using radio waves penetrates the epidermis and heats the dermis. The generated electromagnetic waves and thermal energy change each other, and the generated thermal effect promotes the production of collagen and reduces fine lines.
3. Red light care increases the skin's capacity to penetrate the active ingredients contained in the care product. Daily use can effectively improve skin elasticity and reduce discoloration around the eyes.
4. High-frequency vibration massage supports blood circulation and better absorption of the skin by the active ingredients contained in the care product.

Operating instructions

Power Button – Long press the power button to turn on, short press to switch working mode.

1. RF + red light + vibration (intensity unchanged)
2. RF + red light
3. RF + vibration (intensity unchanged)
4. RF radio frequency

Level switch key: switch three gears (gear only adjusts RF intensity)

Using instructions

Take an appropriate amount of eye cream and spread evenly on the skin care area.

Long press  to turn on device, short press  to switch working mode, and short press  to adjust the RF intensity as needed.

Hold the middle and upper part of the body and gently glide massage probe around eye contours as shown in the picture.

Note: This machine is not waterproof, please do not rinse with water.

Precautions

1. Follow the recommended user manual. If your eyes have inflammation such as conjunctivitis, you should not use this device until it has healed completely. Do not use this product while wearing contact lenses.
2. Remove all metal parts from the facial area before using the massager, such as earrings, necklaces.
3. Do not use the B-Looky PRO while it is charging.
4. Do not repair, disassemble or modify the device on your own. Any irregularities should be assessed by the service.
5. Please do not drop or hit the device to avoid malfunction.
6. When the power cord is damaged, please do not use it to avoid short circuit or fire.
7. Do not use the device in people:
 - with infectious diseases — with skin diseases
 - with heart diseases — pregnant / breastfeeding
 - with high fever — with excessive sunburn
 - with a malignant tumor
 - with acute inflammation
 - with a pacemaker

Storage

1. Apply eye area care products evenly before each treatment to avoid excessive accumulation of eye area products in the massage head.
2. After completing each treatment, wipe the product with a dry towel.
3. Do not wash the device under water or submerge it. The device is not waterproof.
4. Do not use corrosive agents to clean the device.
5. Store it in a cool, dry place.
6. Do not use the device if it is damaged.
7. Do not press the device with great force against the skin.

Charging method

1. When it has low power, please turn off the instrument and connect the charger.
2. During the charging process, the battery symbol on the display flashes, and the battery symbol stops flashing when charging is completed.
3. Please charge it occasionally to avoid aging if you do not use it for long time.
4. Please unplug the charger after charging is finished.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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Beauty Factor Sp. z o.o.
Plac Bankowy 2, 00-095 Warszawa

English

Hereby, Beauty Factor Sp. z o.o. declares that the device equipment Beautifly B-Looky PRO is in compliance with Directive 2014/30/UE.

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beauty Factor Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
4. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.

In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address serwis@beautifly.eu

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I represent that I have read and accepted the
terms and conditions specified in the Warranty
Certificate.